

Revolution Plus Sample Prep System

Firmware Update Procedure – v1.0.5

IMPORTANT: If you have been provided with a Research Mode USB drive, installing Research Mode will replace Kit Mode on the system. To return to Kit Mode, reinsert the Kit Mode USB drive and reboot the system. The same process is used to switch between modes.

Step 1 Press the “O” side of the power switch on the back of the system to turn the system off.



Step 2 Insert the Kit Mode or Research Mode USB drive for the desired operating mode into the USB port on the lower-right side of the back of the system.



Step 3 Press the “I” side of the power switch to turn the system on.



Step 4 The screen will display the firmware version. Verify that the desired mode is displayed, then press “Update.”




Step 5 Do not remove the USB drive or turn off the system while the update is in progress



Step 6 When the update is complete, the screen will prompt you to remove the USB drive. Remove the USB drive.



Step 7 Confirm that the firmware version appears in the lower-right corner of the display. This indicates that the firmware update is complete.

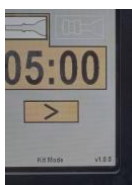

Update Complete

Your Revolution Plus system is now updated and ready for use.

Research Mode users: To switch back to Kit Mode, turn off the system, insert the Kit Mode USB drive, and restart the system.